



RotationVR- Rehabilitation

Module tutorial



Table of contents

1. GENERAL.....	3
1.1. Description	3
1.2. Required accessories	3
1.3. Patient setup	3
2. Settings	3
2.1. Session settings	3
2.1.1. Settings	4
3. Module	4
3.1. Session.....	4
3.2. Shortcuts	5
3.3. Results	7
3.3.1. Report and graphs	7



1. GENERAL

1.1. Description

The **RotationVR-Rehabilitation** software is an assistant for using the rotating chair. This exercise will make it possible to record and store a variety of information during a session.

1.2. Required accessories


The **RotationVR-Rehabilitation** software is an assistant for using the rotating chair. This exercise will make it possible to record and store a variety of information during a session.

1.3. Patient setup

Patient sitting on the rotatory chair.

2. Settings

2.1. Session settings

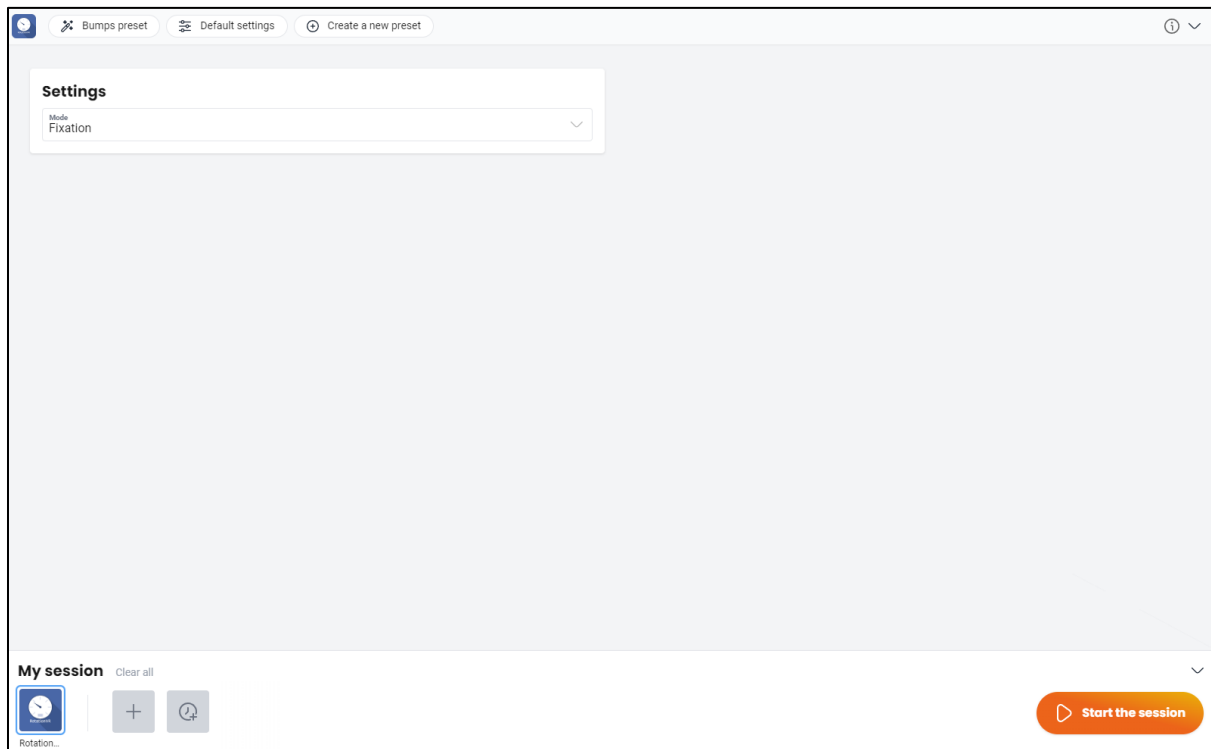

The following setting can be changed before the session:

2.1.1. Settings



Mode:

Selects the mode of use: **Fixation** or **Vection**.

3. Module

3.1. Session



Once the presets have been selected, click on "**Start the session**" in the bottom right corner of the screen.





Each of the rotations performed during the exercise takes place in two stages:

1. The first stage consists of rotating the chair, with the software recording the data relating to the rotation. The user can check the speed of rotation and its progress at any time on the interface.
2. The second phase starts when the rotation is completed, and the chronometer is started. During this time, the user must count the nystagmus and enter it in the related fields. The speed of the slow phase can also be entered in the corresponding field.
3. At the end of the analysis, stop the chronometer before starting the new rotation. When the chronometer is stopped, the fixation or vection time is recorded.

3.2. Shortcuts

During the session, the shortcut list is found by clicking on the Xbox controller icon in the upper right corner of the screen.

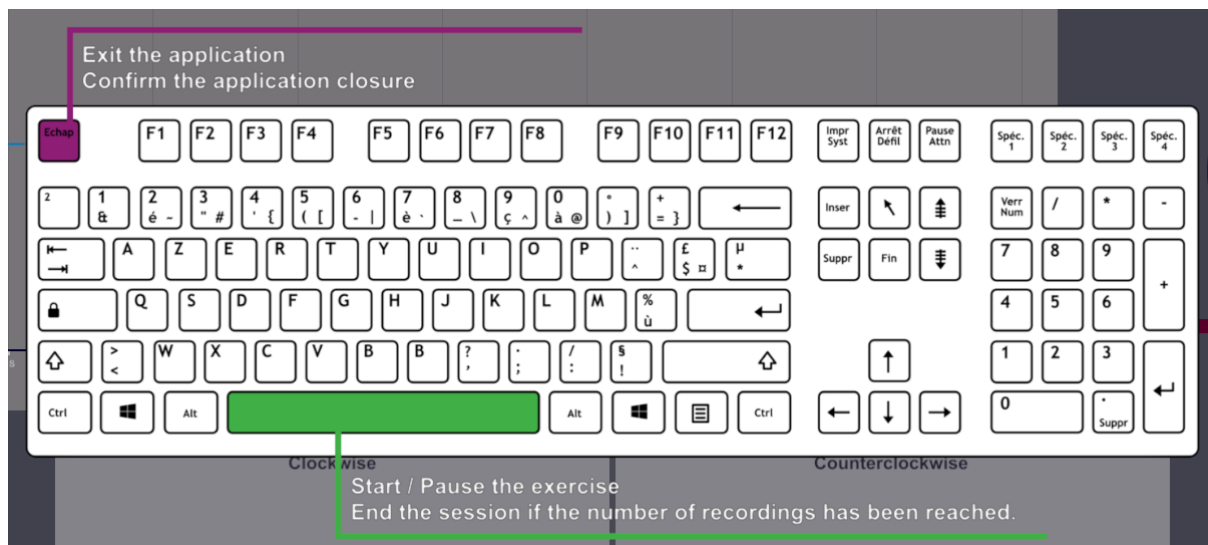

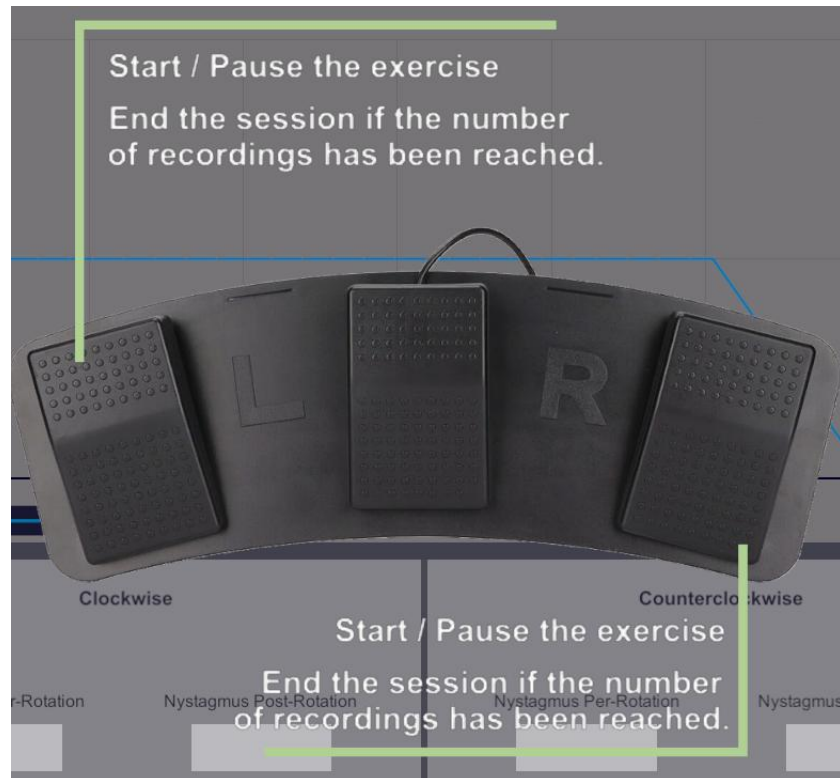

3.3. Results

Once the session is over, you can access the results.

3.3.1. Report and graphs

Click on the histogram icon to access detailed results and the session report.



Session details
22/02/2025 15:58

RotationVR-Rehabilitation

Parameters

Results

Notes

Name	Result			
Rotation direction	1.00	0.00	0.00	1.00
Number of turns	1.00	1.00	1.00	1.00
Time	5.66	5.48	4.78	4.61
Average speed	82.40	66.50	66.50	77.80
VPL	0.00	0.00	0.00	0.00
FrequenceNystagmus	0.00	0.00	0.00	0.00
Minimum acceleration	3449.40	61.20	61.20	86.70